



## Download The Medesto App

Many patients struggle with proper medication adherence. Noncompliance can result in unwanted health complications and contributes to misuse and waste. That's where we come in.

Scan the QR code below for easy access to the Medesto App.



[perigonhealth.org/FAQ](http://perigonhealth.org/FAQ)

844.698.2533

# User Guide

## Medesto Connect

Version C1.0.2



SCAN FOR PDF



## Introduction

Welcome! Thank you for choosing Medesto. Our goal is to empower individuals with health goals or chronic conditions to live better and healthier lives.

**Adherence:** the extent to which patients take medication as prescribed by their doctors.

**Non-Adherence:** the intentional or unwitting failure to take medications as prescribed.

**Many patients struggle with or do not comply to their medication adherence regimen. This non-adherence can have varying impact on a patient's health.**

That's where we come in.

Our innovative ecosystem includes the **Medesto Connect, Medesto Dispenser, Medesto Pods, and the Medesto App**, which work as a connected system to improve adherence and instill peace of mind in patients and caregivers.

The Medesto Connect is a compact design to sit on a countertop, end table, bed side or taken on the go, to provide connectivity and audio and visual reminders. Scan NFC smart tagged medications or cards to indicate usage taken that day. The Medesto Dispenser is designed to dispense solid medications from the pods and when connected to WiFi and the app, tracks the doses taken or missed.

Medesto Pods can be organized and filled at our in-house pharmacy, Perigon Pharmacy 360, and sent directly to patient or local pharmacies with Medesto Pod filling, contact your doctor to transfer your pharmacy.

The Medesto App makes adherence easy with notifications when it's time for a scheduled dose, as well as alerts for missed doses.

## Support and Contact Us

**WARNING:** Contact a healthcare provider before using or changing schedules.


Keep this product and device out of reach from children. Adult supervision required.

Do not exceed suggested use unless directed by a doctor or healthcare.

For any questions regarding the Medesto Dispenser, Medesto App, Medesto Connect, Medesto Pods, Prescriptions or Refills, etc., please contact our Customer Care Team Experts immediately before use.

[perigonhealth.org/FAQ](https://perigonhealth.org/FAQ)  
**844.698.2533**

## Set Up & Use

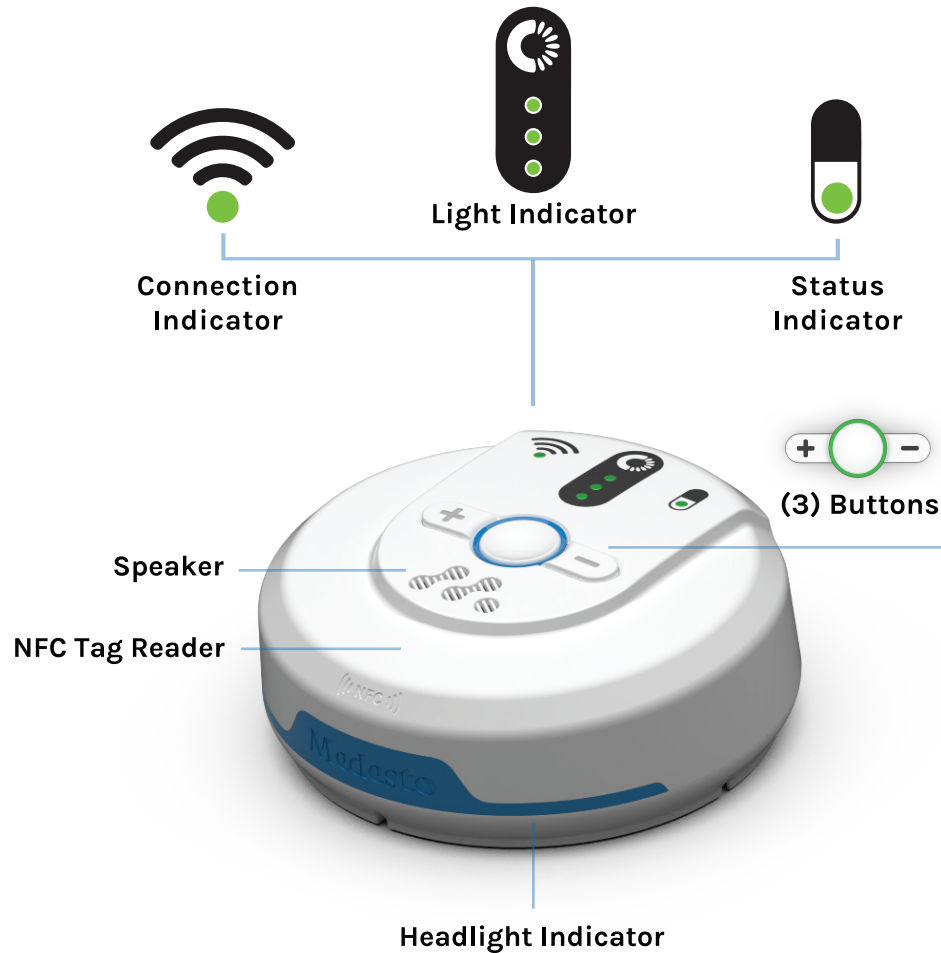

*\*Medesto Connect has no microphone to listen only a speaker to play audio*

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase separation between the equipment and receiver.
- Connect the equipment to an outlet on a different circuit from which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

*Changes or modifications not expressly approved by the party responsible for adherence could void the user's authority to operate the equipment.*

**Medesto™**

MANUFACTURED FOR PERIGON HEALTH 360.  
ASSEMBLED IN CHINA.



US PAT. NO. 10,759,594  
US PAT. NO. 10,279,985  
US PAT. NO. 10,252,843  
US PAT. NO. 11,014,721  
ADDITIONAL PATENTS PENDING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

AUS PAT. NO. 2,019,204,094  
AUS PAT. NO. 2,014,348,856  
EP PAT. NO. 3,068,274  
CA PAT. NO. 2,929,775  
CN PAT. NO. 106,028,887  
JP PAT. NO. 6,598,214  
JP PAT. NO. 6,928,318  
KR PAT. NO. 10-2378899  
SG PAT. NO. 10,201,804,028P



**FCC ID: 2A8H7-LB300-5000**

**CONTAINS FCC ID: 2AC72-ESPWROOM32 AND XMR2020BG95M1  
FCC TESTING WAS CONDUCTED AFTER ASSEMBLY**

**USE ONLY WITH ADAPTER TECH, MODEL: ATSO36T-W240U  
REFER TO USER GUIDE BEFORE USE FOR INFORMATION REGARDING  
SAFE DEVICE USAGE - FOR RESIDENTIAL USE ONLY**

# Medesto Connect User Interface

## Color Indicators



Schedule to be Recorded



Last Schedule Recorded



Last Sch. NOT Recorded/Error



Syncing



Action Needed/Taken



No Action/Sleep/Updating

## Status Indicator



Schedule to be Recorded



Last Schedule Recorded



Last Sch. NOT Recorded



Sleep/ No Status

## Connection Indicator



Connected to Internet WiFi



Connected to Cellular LTE



NOT Connected

## NFC Reader

Place NFC Stickered items; Medications, Pill Containers, or Smart Cards on top of Medesto Connect shown here on the blue area. This will log a current status/scheduled event.

Wait for the lights and sounds to indicate a successful scan.



# Warranty / Disclaimer / FCC Information

THIS LIMITED WARRANTY COVERS ONLY DEFECTS IN WORKMANSHIP AND MATERIALS FOR A PERIOD OF ONE (1) YEAR FROM THE DATE OF PRODUCT PURCHASE AS EVIDENCED BY THE ORIGINAL PURCHASE RECEIPT (THE "WARRANTY PERIOD"). DURING THE WARRANTY PERIOD, PERIGON HEALTH 360, WILL REPLACE OR REPAIR ANY PRODUCT, IN ITS SOLE DISCRETION, THAT IS DEFECTIVE AS A RESULT OF DEFECTS IN MATERIAL OR WORKMANSHIP, AND WHICH IS RETURNED TO PERIGON HEALTH 360 OR TO ONE OF ITS EXPRESSLY AUTHORIZED PARTNERS OR AGENTS BY THE PRODUCT'S ORIGINAL PURCHASER. THE FOREGOING REMEDY IS THE SOLE AND EXCLUSIVE REMEDY FOR ANY BREACH OF WARRANTY AND IS AVAILABLE SOLELY TO THE ORIGINAL PURCHASER OF THE PRODUCT. IN ALL CIRCUMSTANCES PERIGON HEALTH 360'S LIABILITY IS LIMITED TO REPAIR OR REPLACEMENT OF THE PRODUCT THAT ARE RETURNED TO PERIGON HEALTH 360 WITHIN THE WARRANTY PERIOD. IN NO EVENT SHALL PERIGON HEALTH 360'S LIABILITY EXCEED THE VALUE OF THE PRODUCT SOLD. TO ACCESS THIS WARRANTY PLEASE CONTACT PERIGON HEALTH 360 DIRECTLY AT THE PHONE NUMBER PRINTED WITHIN THIS USER GUIDE OR ON YOUR PRODUCT'S PACKAGING.

THIS WARRANTY DOES NOT APPLY TO ANY DAMAGE, MALFUNCTION, DEFECT, LOSS, HARM OR OTHER CIRCUMSTANCE CAUSED BY USE, ABUSE, ACCIDENT, MISUSE, ELECTRICAL OR POWER FAILURE, INABILITY TO USE THIS PRODUCT, OR DISSATISFACTION WITH THE PERFORMANCE OF THE PRODUCT. IT ALSO DOES NOT COVER ANY USE THAT IS NOT IN ACCORDANCE WITH THE USER GUIDE AND MANUFACTURER SPECIFICATIONS, NOR ANY COMMERCIAL USE OF THE PRODUCT. THIS WARRANTY SPECIFICALLY EXCLUDES COVERAGE OF PRODUCTS FOR WHICH PERIGON HEALTH 360 HAS NOT RECEIVED PAYMENT OR WHICH ARE USED BY OR IN THE POSSESSION OF ANY PARTY OTHER THAN THE ORIGINAL PURCHASER.

THIS WARRANTY GIVES THE ORIGINAL PURCHASER OF THE PRODUCT SPECIFIC LEGAL RIGHTS WHICH MAY VARY FROM STATE TO STATE (PROVINCE TO PROVINCE IN CANADA). PERIGON HEALTH 360'S RESPONSIBILITY FOR DEFECTS IN MATERIAL AND WORKMANSHIP OF THE PRODUCT SHALL BE LIMITED TO REPAIR AND/ OR REPLACEMENT OF THE PRODUCT AS SET FORTH IN THIS WARRANTY. NOTWITHSTANDING ANY TERM HEREIN, ANY AND ALL WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY WARRANTY OF MERCHANTABILITY OR WARRANTY OF FITNESS FOR A PARTICULAR PURPOSE ARE HEREBY EXPRESSLY DISCLAIMED.

TO THE MAXIMUM EXTENT PERMITTED BY APPLICABLE LAW, IN NO EVENT SHALL PERIGON HEALTH 360, INCLUDING ITS SUBSIDIARIES, AFFILIATES, PARTNERS, REPRESENTATIVES, AND AGENTS BE LIABLE FOR ANY INCIDENTAL, SPECIAL, INDIRECT OR CONSEQUENTIAL DAMAGES, WHETHER RESULTING FROM THE USE, MISUSE, OR INABILITY TO USE THIS PRODUCT, FROM DEFECTS IN THE PRODUCT, OR RELATING TO OR ARISING OUT OF ANY CLAIM, DEMAND, LIABILITY, CAUSE OF ACTION OR OTHER LEGAL ACTION, INCLUDING THOSE ALLEGING PERSONAL INJURY, DEATH, LOSS OF PROPERTY, OR ANY OTHER CLAIM WHETHER IN CONTRACT OR IN EQUITY.

PLEASE NOTE THAT SOME STATES (PROVINCES IN CANADA) MAY NOT ALLOW THE EXCLUSION OF INCIDENTAL, SPECIAL, INDIRECT OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATIONS MAY NOT APPLY TO YOU.

THE PRODUCT IS OFFERED IN CONDITION UPON THE PURCHASER'S AND/OR USER'S ACCEPTANCE OF ALL TERMS, CONDITIONS, AND NOTICES WITHOUT MODIFICATION OF THE TERMS, CONDITIONS AND NOTICES CONTAINED HEREIN. MODIFICATION OF THE PRODUCT THROUGH USE OR MISUSE, APPLICATION OF ADDITIONAL COMPONENTS AND/OR REMOVAL OF PROVIDED COMPONENTS, CAN LEAD TO PRODUCT MALFUNCTION CAUSING SERIOUS RISK. RISKS INCLUDE BUT ARE NOT LIMITED TO PRODUCT DAMAGE, SERIOUS PHYSICAL INJURY, DISABILITY AND/OR DEATH. PURCHASER AND/OR USER ARE RESPONSIBLE FOR ANY MODIFICATIONS PERFORMED ON THE PRODUCT BY PURCHASER, USER, AND/OR ANY THIRD PARTY.

PERIGON HEALTH 360 IS NOT RESPONSIBLE AND DISCLAIMS ALL LIABILITY AND RESPONSIBILITY FOR ANY MODIFICATIONS DONE TO THE PRODUCT THROUGH THE USE OR MISUSE OF THE PRODUCT, INCLUDING ADDING OR REMOVING COMPONENTS THAT ALTER THE FUNCTIONALITY OF THE PRODUCT AND/OR IMPROPER INSTALLATION OR ADJUSTMENT OF THE PRODUCT. THIS PRODUCT MAY NOT BE USED FOR UNLAWFUL PURPOSES AND SUCH USES ARE EXPRESSLY PROHIBITED UNDER THE TERMS AND CONDITIONS OF ITS USE.

IF YOU DO NOT ACCEPT THESE TERMS, CONTACT US WITHIN 14 DAYS OF PURCHASE. PERIGON HEALTH 360 RESERVES THE RIGHT TO CHANGE, MODIFY, OR AMEND THE TERMS, CONDITIONS AND NOTICES UNDER WHICH ITS PRODUCTS ARE OFFERED AT ANY TIME WITHOUT NOTICE.

## Home Button



Schedule to be Recorded



Last Schedule Recorded



Review Last Schedule/ Error



Syncing



Action Needed



Updating/ Sleep Mode\*

\*Awake device from the power saving, Sleep Mode, by holding the home button for **3 seconds**


## Light Indicators



Battery 100%



Battery <70%

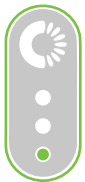

Battery <30%



Battery <10%



Use Buttons



Pressed



Pressed

## Arrow Buttons

1. Plus Button



Yes  
True  
Well  
Correct  
Always  
At Least Once  
Option A  
Volume Up  
Brightness Up



2. Minus Button

No  
False  
Not Well  
Wrong  
Never  
None  
Option B  
Volume Down  
Brightness Down

### Listen Carefully

Take Action, press the purple center button to listen to audio messages and answer audio questions.

## Tech Specs

### Functional Specifications

	24DVC Use Only with AC Adapter
Input Voltage	Adapter Tech, Model: ATSO36T-W240U  AC Power Adapter to be used only plugged into a standard power outlet. (Overvoltage Category II)
Input Power	1.3A Use Only with AC Adapter Adapter Tech, Model: ATSO36T-W240U
Input Frequency	DC 50-60Hz Use Only with AC Adapter Adapter Tech, Model: ATSO36T-W240U
Power Fail	Battery Back up
WiFi	802.11 b/g 2.4Ghz
Cellular SIM	802.15 LTE Cat M1/Cat NB1/EGPRS
Duty Cycle	Continuous Use
EMC Compliance	FCC Part 15, Class B
Safety	IEC 61010-1
Data Security	HIPAA



 **Medesto**